

Instruction of using PyXAS GUI

Reference: “PyXAS-an open-source package for 2D Xray near-edge spectroscopy analysis”, J. Synchrotron Rad. 2020, 27, 567

Source code: <https://github.com/gmysage/pyxas>

Revised on 2020/6/4, by Mingyuan Ge (mingyuan@bnl.gov)

Contents

Individual XANES image fitting	2
Sample files	2
Step-by-step of 2D XANES fitting.....	2
Data loading (Tab: Prep.) (see Fig. 1):	2
Image alignment (Tab: Prep.) (see Fig. 1):	2
Edge normalization (Tab: Prep. & Norm.) (Fig. 2, Fig. 3)	3
2D XANES fitting (Tab: Fit Spec.) (Fig. 4)	4
Segmentation/Mask (Tab: Img. Tools).....	5
Miscellaneous	5
Colormix (Tab: Img. Tools) (Fig. 5).....	5
Basic image processing (Tab: Img. Tools) (Fig. 5):	6
Fit XANES peak position (Tab: Others) (Fig. 6).....	6
Batch processing (e.g., 3D XANES).....	8
Data preparation	8
Step-by-step	8
Load files	8
TXM normalization	9
Edge-normalization.....	9
Image alignment	9
Peak regulation	9
Fitting	9
Mask	9
Color mix	10
Reference and energy list.....	10
Save parameter and batch fitting	10
File saving.....	10

Edge normalization (Tab: Prep. & Norm.) (Fig. 2, Fig. 3)

1. Take -log of image if necessary. For the current example, we don't need to take negative natural log. For regular TXM XANES image, negative natural log is required, by clicking "[Norm. TXM \(-log\)](#)" button.
2. Check energy range of pre-edge and post-edge (**Fig. 2**)
Select the one of the ROI listed in the "[ROI list](#)", click "[Plot Spec.](#)". For this example (see ROI spectrum below), pre-edge can be defined with energies below 8.33 keV, and post-edge can be defined for all energies above 8.4 keV.

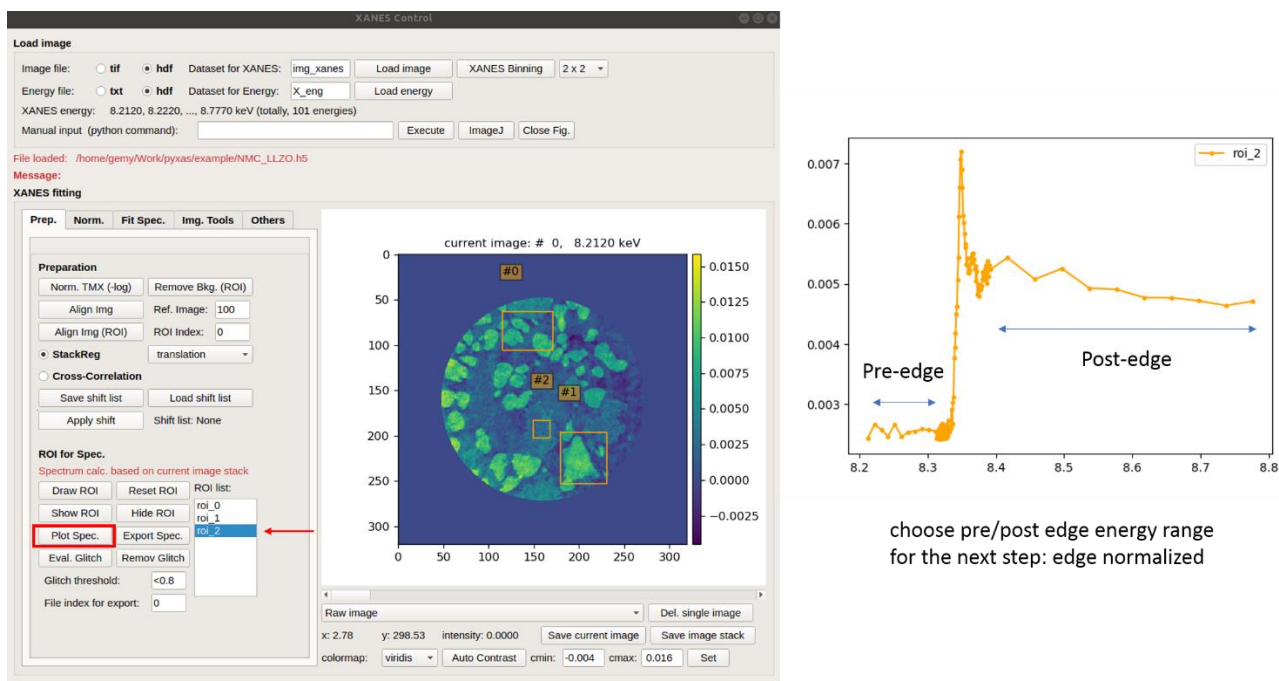


Fig. 2. Draw ROI and plot the spectrum, determine the energy range for pre/post absorption edges.

3. Under the "[Norm.](#)" tab, input the pre and post-edge energy range in the relevant fields: "[Pre-edge start/end](#)" and "[post-edge start/end](#)" (**Fig. 3**)
4. To check how does the edge normalization work, click "[Norm Spec \(ROI\)](#)". It will plot the raw spectrum and edge-normalized spectrum averaged from the ROI region as selected.
Note: if this is material that will be used as a reference material, click the "Save Spec" button to save the normalized spectrum for use as the reference spectrum required for the subsequent XANES fitting. (**Fig. 3**)
5. Once we are satisfied with the defined energy range of pre/post edge, we can perform the edge normalization for the whole image, simply by clicking "[Norm Image](#)" (**Fig. 3**)
Note: we provide two methods for edge normalization.
 - i. If image is noisy, use method 1. It will calculate slope of the absorption curve of pre/post range by averaging all the pixel values from areas that contains materials and keep the slope unchanged and then fit the offset of the absorption curve for individual pixel.
 - ii. If image has good signals, use method 2. It will calculate the linear slope and offset of absorption curve individually for each pixel.

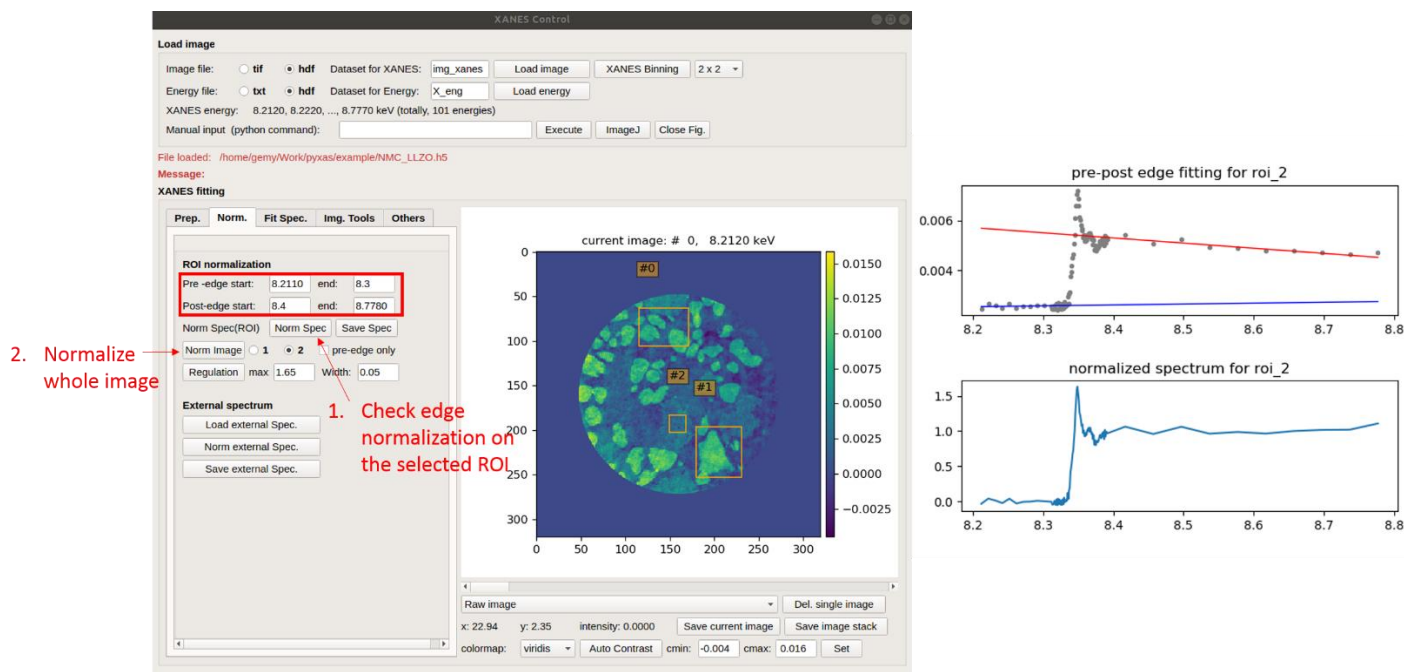


Fig. 3. Edge normalization on selected ROI (e.g., ROI2) and the whole image. Left: screen shot of package interface. Right: edge-normalization for ROI2

2D XANES fitting (Tab: Fit Spec.) (Fig. 4)

1. Load reference spectra: click “Load Ref.” to load “ref_NiO.txt” and “ref_LiNiO2.txt”
We can assign labels for these two ref. spectra. E.g, input “Ni2, Ni3” in the “Elem.” filed.
2. Click “Fit 2D”
Note 1: The fitting will use the image data which is currently displayed in the image display window. Make sure we selected the proper dataset. Mostly, we should use “Image updated”.
Note 2: for fitting using iterative method (e.g, Conj. Grad), give the relevant parameters, and then click “Fit 2D (iter)”
3. After fitting, additional image datasets will be added in the image display window: “XANES Fit thickness”, “XANES Fit (ratio, summed to 1)”, “XANES Fit (Elem, concentration)” and “XANES Fit error”.
Note: $\text{XANES Fit (Elem, concentration)} = \text{XANES Fit (ratio, summed to 1)} \times \text{XANES Fit thickness}$.
4. Fitting evaluation
We can view the “XANES Fit error” to check the fitting quality.
We can also check the fitting results from specific ROI region. E.g., if we would like to check ROI 1, put “1” in the “ROI #” field. Then click “Plot ROI fit”. It will plot both the experimental data and fitted curve. A good match indicates a good fitting.
5. The fitting results can be saved in a .h5 file by clicking “Save 2D Fit”.

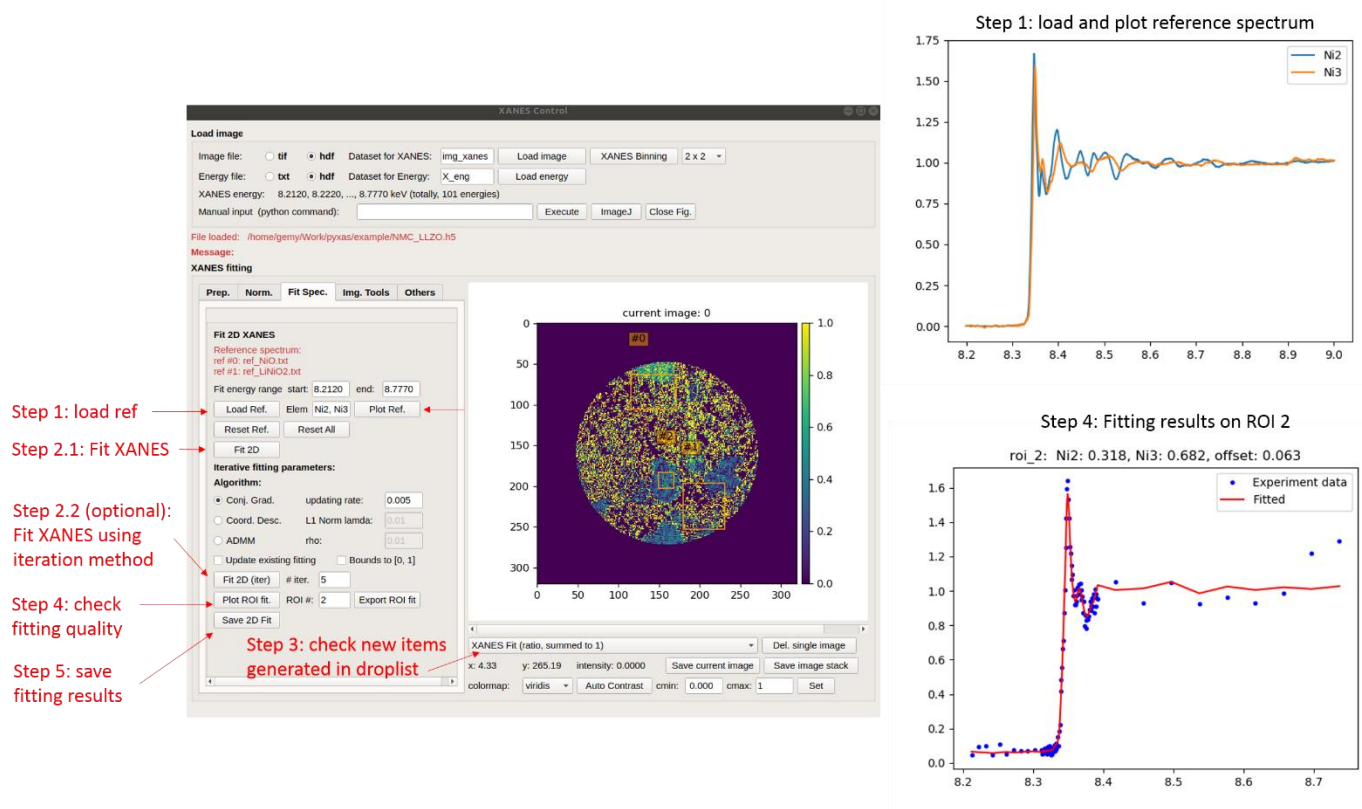


Fig. 4. XANES fitting

Segmentation/Mask (Tab: Img. Tools)

1. Threshold mask:

E.g., generate mask from “XANES Fit error”: display the “**XANES Fit error**” and select a threshold value. E.g., we would like to remove regions with fitting error larger than 0.3. Then we can write “>0.3” in the field of “**Gen. Mask1(2) Thresh:**”, then click “**Gen. Mask1(2)**” button.

Note: the mask will be applied to all the images in the package.

To remove the mask, simply by clicking “**Rmv, Mask1(2)**”

2. Clustering mask:

- Choose the “Raw image” in the image display window.
- Under “**Clustering Mask**”, check the checkbox “**Use stack**”. If we would like to use part of image stack, e.g., from image #20 to #40, put “20” in “**start:**” and “40” in “**end:**”.
- Put the desired number of components, e.g., “4” in the field of “**comp #**”
- Click “**Gen. Mask**”
- A new image dataset “Smart Mask” will be added to the image display window for visualization.
- For the 4 masks generated (because we set comp# to be 4), if we would like to apply the second mask to the image, in the image display window, navigate the slider-bar to the second mask, then click “**Apply mask**”. It will propagate the mask to all images. We can remove it by click “**Rmv. mask**”.

Miscellaneous

Colormix (Tab: Img. Tools) (Fig. 5)

After XANES fitting, we can use colormix to generate a colored image to represent the distribution of element with different oxidation states, e.g., Ni^{2+} and Ni^{3+} .

- In the droplist in image display window, select the fitted results, e.g., “**XANES Fit (Elem, concentration)**” or “**XANES Fit (ratio, summed to 1)**”
- Under section of “**Img. Tools**”, click “**Color Mix**”. The colored image will be added to image display window. We can save the color image by clicking “**Save current image**” to save it.

Basic image processing (Tab: Img. Tools) (Fig. 5):

We provide basic image processing tools, such as [image dilation](#), [erosion](#), [fill-hole](#), [median filters](#) to remove noises.

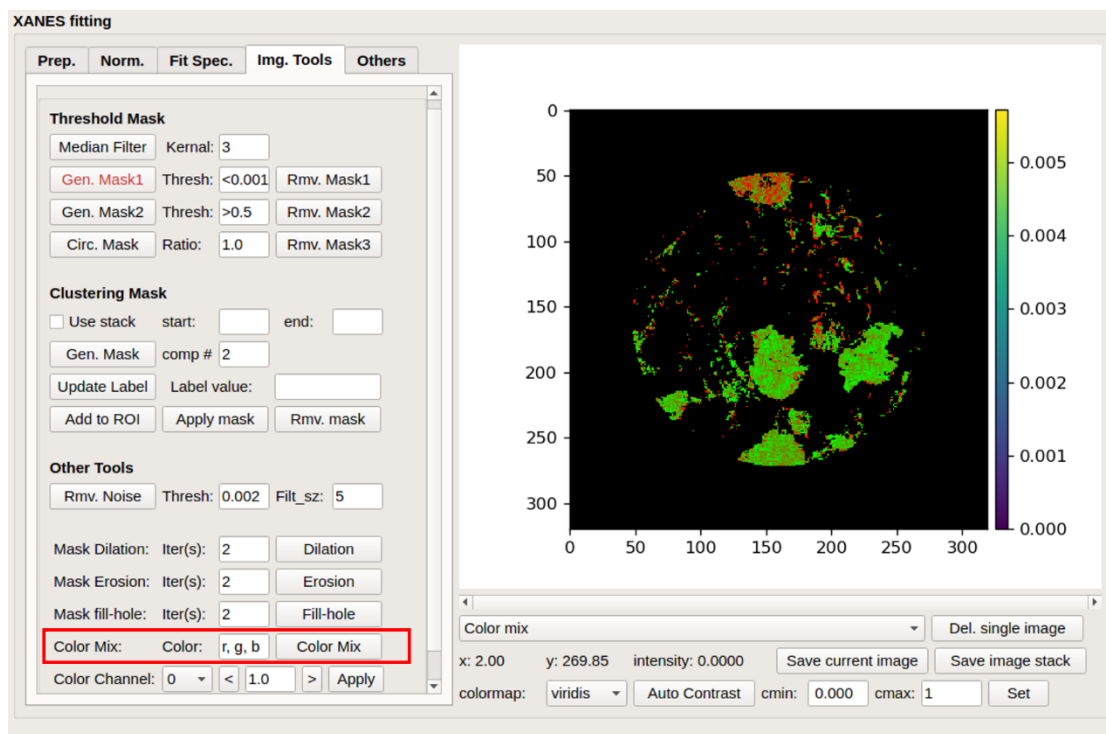


Fig. 5. Image tools

Fit XANES peak position (Tab: Others) (Fig. 6)

In case we don't have reference spectrum at hand, it is convenient to look at the peak position of XANES spectrum as a first approximation to the variation of chemical state (oxidation state).

1. In the droplist in image display window, select "Image update"
2. Under section of "Other", input the estimated peak position range. See Fig. 6 left.
3. Here, it uses spline curve to fit individual spectrum. It is recommended to test the fitting parameter "smooth" and "order" by fitting ROI spectrum first.
 - a. Input the ROI # (assume you have created ROI already, check Fig. 1)
 - b. Click "Fit ROI", it will generate a figure showing the fitting results and indicate the peak position. (Fig. 6 right)
4. To fit the whole image, click "Fit image"
5. Fitted image will be added to the droplist in the image display window
6. (optional) it can scale the fitted image to (0, 1) by assigning the min/max value of peak position.

Note: if the "Image update" has very low intensity, e.g., pixel value around 0.001, it is recommended to scale the image intensity by e.g, x1000.

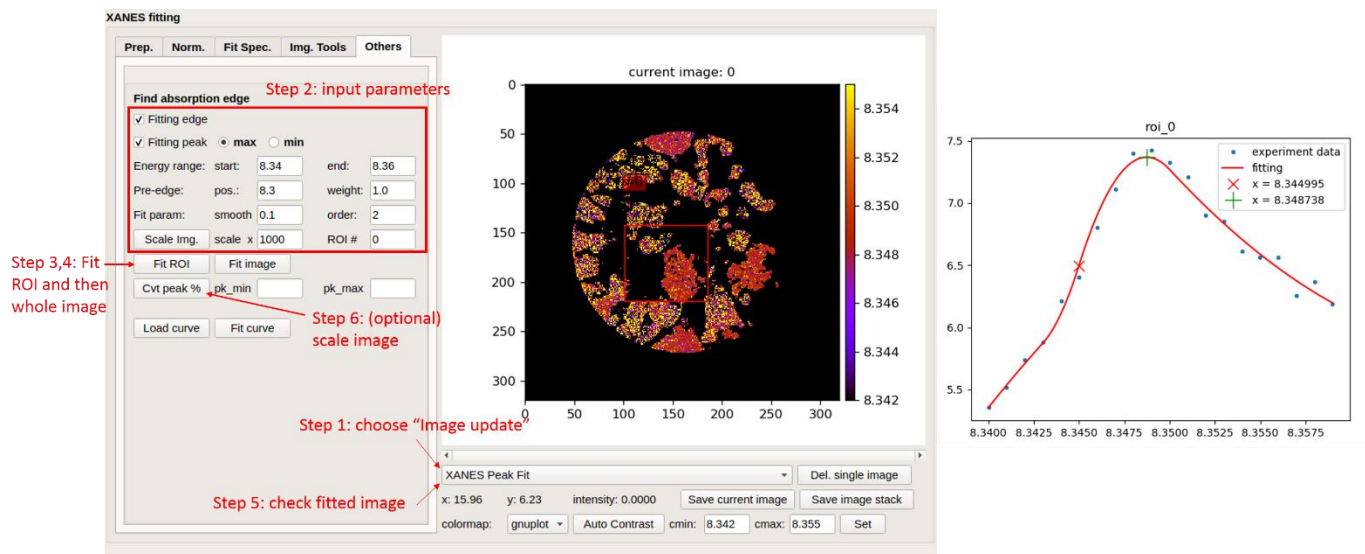


Fig. 6. Fitting peak position of XANES spectrum

Batch processing (e.g., 3D XANES)

Lunch program: [python batch_processing.py](#)

The 'Fitting Parameter' window is a graphical user interface for configuring XANES data fitting. It includes sections for file selection, parallel computing, normalization, peak regulation, image alignment, fitting parameters, cluster mask comparison, color mix, and reference spectrum loading. The 'Files' section shows a file path and a message indicating 540 files were loaded. The 'Parallel computing' section shows 1 CPU. The 'Edge normalization' section shows 'Norm absorp. edge' set to 'Yes' and 'Global slope' set to 'No'. The 'Image alignment' section shows 'Image alignment' set to 'No' and 'Reference index' set to '-1'. The 'Fitting' section shows 'Fit pre-edge' and 'Fit post-edge' checked, with 'Eng. range' from 8.2 to 8.7. The 'Cluster Mask comp.' section shows 'Cluster Mask comp.' set to 0. The 'Color Mix' section shows 'Color Mix' set to 'r, g, b'. The 'Reference spec. & X-ray energy' section shows 'Load Ref.' and 'Load XEng' buttons.

Data preparation

Image file:

All 2D xanes image stacks (e.g., xanes1.tiff, xanes2.tiff, ... xanes_1001.tiff) are saved in one folder.

Energy list:

Energy list save in a txt file. Assume all xanes images are collected using same x-ray energy list

Reference spectrum:

2-column data (energy vs. intensity) saved in txt files, e.g., Ref_1.txt, Ref_2.txt.

Step-by-step


Load files

- 1.1. Click “Open” to direct into the file directory, select one of the xanes image stack (e.g., xanes_1.tiff), then click OK.
- 1.2. Based on file type, choose “tiff” or “hdf” accordingly. If it is hdf file, write the attribute name pointing to the image dataset, e.g., “img_xanes”, to the textbox on right side of “hdf”
- 1.3. Click “Load”. All files loaded with success will be displayed on the “Message” line

This screenshot shows the 'Files' section of the 'Fitting Parameter' window. The 'File path' is '0Q1/NMC811/20cy/xanes_assemble' and the 'Open' button is visible. The 'File name start with' is 'xan', 'format' is '.tiff', and 'Load' button is visible. The 'File type' is set to 'hdf' and 'img_xanes' is entered in the adjacent text box. A message at the bottom states: 'Message: 540 files loaded: [xanes_2D_slice_000.tiff ... xanes_2D_slice_539.tiff]'.

TXM normalization

Depending on type of data to choose whether need to take $-\log()$. E.g., for regular TXM xanes, we need to take $-\log()$. For 3D XANES tomography (e.g., 3D reconstructed data taken at different x-ray energies) does **NOT** need to take $-\log()$.

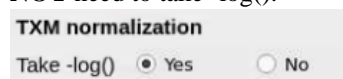


TXM normalization

Take $-\log()$ ☒ Yes ☐ No

Edge-normalization

Choose “yes” if need to perform edge normalization. Input the energy range for pre-edge and post-edge. Unit is keV.

Note: If choose “Global slope”, it is equivalent to use edge normalization method “1” in the previous related section in 2D xanes analysis using xanes_2D.py. (page 3)

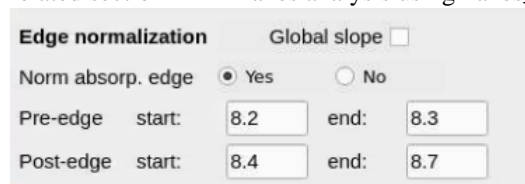

Edge normalization ☐ Global slope ☐

Norm absorp. edge ☒ Yes ☐ No

Pre-edge start: end:

Post-edge start: end:

Image alignment

Choose “Yes” if need to align the image stack. Input the which image will be used as reference. “-1” means the last image in each image stack will be used as reference, and all other images in the image stack will be aligned with respect to it.

You can also apply a circular mask to the images, using parameter “Ratio”. Ratio=0.9 means a circular mask with diameter of 0.9x image width will be applied


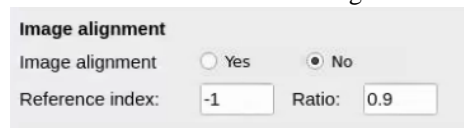


Image alignment ☐ Yes ☒ No

Reference index: Ratio:

Peak regulation

Most of time, we don't need it.

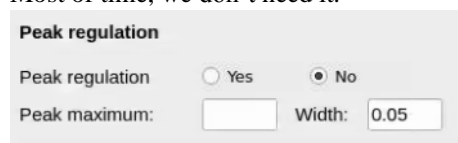

Peak regulation ☐ Yes ☒ No

Peak maximum: Width:

Fitting

Parameters used for fitting. If “Fit iter” is checked, the iteration method will be applied after regular xanes fitting. The parameter “lambda” controls the smoothness of fitting results. If image signal is very good, set lambda=0.



Fitting

☒ Fit pre-edge ☒ Fit post-edge

Eng. range start: end:

☒ Fit iter

Updating rate: # Iter:

Bounds to: low: high:

lambda:

Mask

1. If “Cluster Mask comp.” > 0 (e.g., 2), it will use clustering method to create mask to segment the image into (e.g., 2) masks. Check the paper ([J. Synchrotron Rad. 2020, 27, 567](#)) for details. The clustering segmentation will be applied based on the aligned images (if image alignment is chosen) otherwise on the very raw images loaded.
2. Threshold mask: mask be calculated by the threshold value on the “thickness image” and “fitting_error” image. See previous “Segmentation” section on page 4.

Cluster Mask comp.

Cluster Mask comp.

Threshold: thickness: error:

Color mix

Two types of colormix images will be generated based on “xanes_fit_ratio” and “xanes_fit_concentration”

Color Mix

Color Mix choose from 'r, g, b, c, p, y'

Reference and energy list

Load the reference spectrum and energy list.

Note: all xanes files must have the same energy list during experiment.

Reference spec. & X-ray energy

Reference spectrum:
 ref #0: .../media/mingyuan/Seagate Backup Plus Drive/TXM_2019/r
 ref #1: .../media/mingyuan/Seagate Backup Plus Drive/TXM_2019/r
 Energy: 101 energies: 8.191998738300539, 8.201999949077775,

Save parameter and batch fitting

MUST save the parameter first, before click “Batch fitting”

You can also load previously saved parameter file, and do “Batch fitting”

File saving

Finally, all fitting files are saved in the subfolder “Fitted”.

Note: keep eye on the output in the command terminal. If you want to terminate the program, using “ctrl + C”